## Your Breath, Our Care.



#### Resvent Medical Technology Co., Ltd.

Add:BC601, BC602, Gaoxinqi Factory, District 67, Xingdong Community, Xin'an Street, Bao'an District, 518100 Shenzhen, PEOPLE'S REPUBLIC OF CHINA Tel: 0086-755-23027370-808 Mail:resvent@resvent.com Https: www.resvent.com



## iBreeze+ Series

Redefine your premium CPAP



# Redefine premium CPAP


#### Innovation that excites

iBreeze+ premium CPAP package is the only high-performance all-around PAP solution in the world.

The relentless pursuit of perfection & care has brought quality & portability together.

#### Innovation that performs


iBreeze+ is an unique circular CPAP/APAP device designed to perform with optimum quality therapy, for the most comprehensive sleep and respiratory disorders. iBreeze+ family introduces the first, CPAP, APAP with the intelligent REG algorithm and detailed reporting platform.

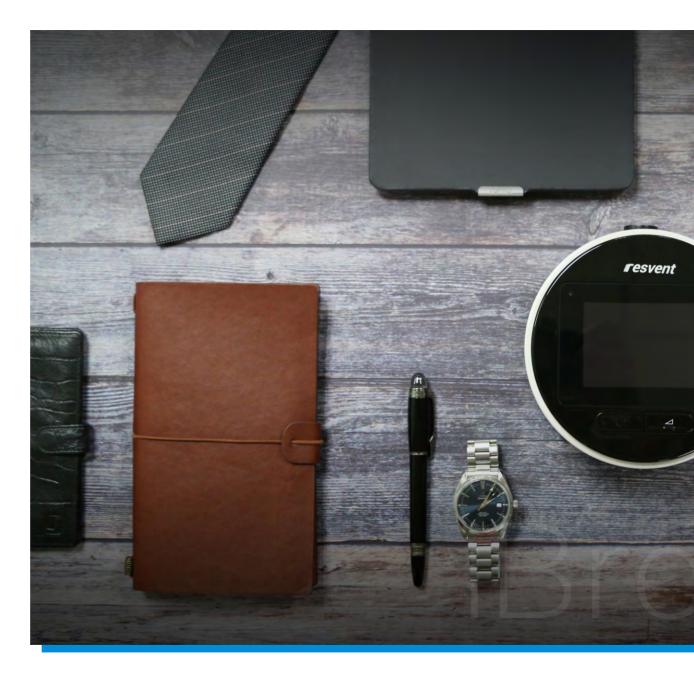



## Portable Comfortable Intuitive Performing Connectivity



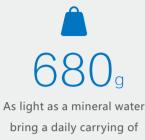

iBreeze+ no any headache

at all.



Same size as your smartphone keeps the consistent decent of your carry bag as before.



travelling.



## Portable

#### Your routine headache, we perceive.

As a regular CPAP user, taking your old therapy device during travelling days may consume you lots of extra attention & energy. A heavy weight, a hulking carry bag and AC dependent system have been depriving the elegance from the days.

#### Your daily convenience, we care.

A portable CPAP as iBreeze+ will meet all your imagination & expectation of a sleep on your go. Without compromise on performance within 680g, iBreeze+ is indeed representing the highest state of art in field.



\*Therapy pressure @ 10cmH<sub>2</sub>O and standard testing condition.





## Comfortable

#### Pleasing at home, and on your go.

Comfortableness means a lot to a regular user of CPAP. Integrated with resvent brilliant experiences in respiration & ventilation, abundant advanced technology and algorithm on pressure release, synchronization and event diagnosis are adopted on iBreeze+ to satisfy even the most harsh customers. A noise level down to 28dB will maintain a quiet environment for all family members throughout the night.

#### Heating tube

iBreeze+ has been designed with integrated power supply system on humidifier and can be connected with heating tube, which will eliminate the water condensation even in most harsh environment.

#### Silencer

With the help of silencer connected with main unit, a quiet and peaceful therapy night becomes reachable for you and families.





#### Intelligent humidifier

iBreeze+ series intelligent humidifier is designed with cover and seal water tank. System will continuously monitor the water level, and stop heating & give a visible alert, while the water level comes low, in avoiding patient's discomfort leaded by dry air and delivery the maximum of safety.




## Intuitive


#### Touch, ease your every setup.

As a daily partner during sleep, easy to use without taking your extra hassle is goal we are chasing at all time. Intuitive touching screen supported by flat UI is available on iBreeze+ series which will bring you a real interaction experiences while using. See, touch, and get it done becomes to possible on a compact system.

#### Ergonomic & thoughtful design

The ergonomic design are everywhere on iBreeze+ series. Braille type of button design, rotatable humidifier outlet connector, and one hand operable humidifier and battery release lock make your daily use with ease and pleasance.

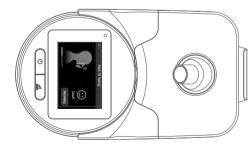

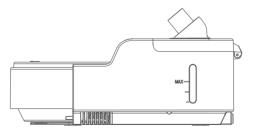




#### Summary on main unit

By accessing our iMatrix cloud system, you could of full scenario of your therapy report. However, a sur daily report still will be generated for your conveni glance every morning.



Clinician's view: informative & supportive


		Information	
		2018-01-01	
	<b>5.6</b> h	<b>6.8</b> /h	
	98%	60 topm	
btain a			
obtain a mmary ient			

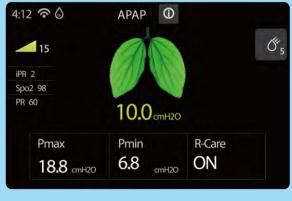

Patient's view: intuitive & clear

## Performing

Compact body, with cutting edge technology.

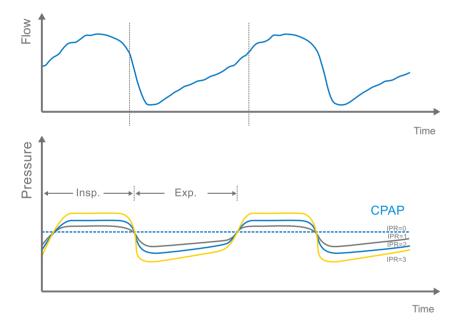

#### IPR

Unique IPR (Intelligence Pressure Release) algorithm during inspiration & expiration offer patients more comfortable experience during ventilation therapy. Three IPR grades are available to meet the personalized requirement.



#### REG

Algorithm innovated by resvent can continuously determine and monitor respiratory events (AHI, OAI, CAI, HI, SNI, FL, RERA, PB, etc), analyze the physical status of patients, and guarantee the delivery of the optimality of therapy pressure based on that at all time.



## Connectivity

#### Access cloud with patient-centric in mind

With integrated wireless module, iBreeze+ can easily link to ResAssist cloud system, a central data gateway. Connect with your smart device and an overall perspective of therapy data can be accessed and viewed easily. A wireless SpO<sub>2</sub> monitoring protocol has been implemented on iBreeze+ to ensure an optimized user experiences.



# ResAs







### ResAssist<sup>™</sup>



Connect with smart phone / tablets will maximize the convenience. You can obtain the overall perspective of patients' data from your PC. Web end offer the most flexible scheme for data management.



## Specification

	Туре	СРАР	APAP
	Model	20C Plus	20A Plus
Modes	СРАР	•	•
	APAP	-	٠
	APAPW		•
	Screen size/inch	3.5	3.5
	Touch screen	•	
	Dimension (w/o humidifier)/mm	270*152*140	
Basic information	Dimension (w/ humidifier)/mm	Ф150*Н68	
	Weight (w/o humidifier)/g	Weight (w/o humidifier)/g 1620	
	Weight (w/ humidifier)/g	680	
	Pressure range/cmH <sub>2</sub> O	4-20	4-20
	Parameter monitoring	•	•
	Event detection	•	•
	Auto Ramp	•	•
	IPR	•	•
	Altitude compensation	•	•
	Leakage compensation	٠	•
	Circuit compensation	•	•
Value added function	Mask fit	٠	•
	Auto illumination	0	0
	Auto Start	٠	•
	Auto Stop	•	•
	Water level monitoring	•	•
	Auto humidification	•	•
	Pre-heating	•	•
	TF card	•	•
Data management	Cloud data service	0	0
Data management	4G modem	0	0
	Wi-Fi & Wireless module	0	0

Standard

Optional

Not available

## Ordering information

Optional accessories	Description	Part number
	Tube (15mm)	50101-010002-00
Tube	Tube (19mm)	50101-010001-00
	Heating tube (Integrated)	50102-000001-00
	Nasal mask	50103-010002-00
Mask	Nasal pillow	50103-000004-00
	Full face mask	50103-010013-00
	Humidifier	60201-000031-00
Humidification	Water chamber	10001-000105-00
	Water less humidifier (HME)	60101-000003-00
Silencer	Silencer (w/o humidifier)	60201-000048-00
Filter	Air filter	10006-000025-00
iPOM SpO <sub>2</sub> module	wireless SpO <sub>2</sub> watch	50302-000003-00
Connection module	Wi-Fi & Wireless module	09901-010003-00
Connection module	Cellular module	09901-010004-00
ExtPowerPack	ExtPowerPack (6700mAh) / 13h battery backup	60201-000077-00
	36 Watt, 100-240VAC, 50/60Hz (w/o humidifier)	05205-000003-00
Power adaptor	65 Watt, 100-240VAC, 50/60Hz(w/ & w/o humidifier)	05205-000002-00
	Vehicle charger	05205-000008-00
	Hard case (for main unit)	70104-000001-00
Carrying solution	Portable accessories case	70104-000004-00
	Comprehensive travelling bag	70103-000004-00



